




Ordering Information

Shape Code	Usable Length	Catalog No. 6Fr	Catalog No. 7Fr	Catalog No. 8Fr	Tip Shape	
JL3.5	100cm	HW60JL350P10000	HW70JL350P10000	HW80JL350P10000	Judkins Left	
JL3.5 ST	100cm	HW60JL350P10001	HW70JL350P10001	HW80JL350P10001		
JL3.5 SH	100cm	HW60JL352P10000	HW70JL352P10000	HW80JL352P10000		
JL3.5 STSH	100cm	HW60JL352P10001	HW70JL352P10001	HW80JL352P10001		
JL4.0	100cm	HW60JL400P10000	HW70JL400P10000	HW80JL400P10000		
JL4.0 ST	100cm	HW60JL400P10001	HW70JL400P10001	HW80JL400P10001		
JL4.0 SH	100cm	HW60JL402P10000	HW70JL402P10000	HW80JL402P10000		
JL4.0 STSH	100cm	HW60JL402P10001	HW70JL402P10001	HW80JL402P10001		
JL4.5	100cm	HW60JL450P10000	HW70JL450P10000	HW80JL450P10000		
JL4.5 ST	100cm	HW60JL450P10001	HW70JL450P10001	HW80JL450P10001		
JL4.5 SH	100cm	HW60JL452P10000	HW70JL452P10000	HW80JL452P10000		
JL4.5 STSH	100cm	HW60JL452P10001	HW70JL452P10001	HW80JL452P10001		
JR3.5	100cm	HW60JR350P10000	HW70JR350P10000	HW80JR350P10000		Judkins Right
JR3.5 SH	100cm	HW60JR352P10000	HW70JR352P10000	HW80JR352P10000		
JR4.0	100cm	HW60JR400P10000	HW70JR400P10000	HW80JR400P10000		
JR4.0 SH	100cm	HW60JR402P10000	HW70JR402P10000	HW80JR402P10000		
AL0.75	100cm	HW60AL070P10000	HW70AL070P10000	HW80AL070P10000	Amplatz Left	
AL0.75 SH	100cm	HW60AL072P10000	HW70AL072P10000	HW80AL072P10000		
AL1.0	100cm	HW60AL100P10000	HW70AL100P10000	HW80AL100P10000		
AL1.0 SH	100cm	HW60AL102P10000	HW70AL102P10000	HW80AL102P10000		
AL1.5	100cm	HW60AL150P10000	HW70AL150P10000	HW80AL150P10000		
AL1.5 SH	100cm	HW60AL152P10000	HW70AL152P10000	HW80AL152P10000		
AL2.0	100cm	HW60AL200P10000	HW70AL200P10000	HW80AL200P10000		
AL2.0 SH	100cm	HW60AL202P10000	HW70AL202P10000	HW80AL202P10000		
SAL0.75	100cm	HW60AL070P10002	HW70AL070P10002	HW80AL070P10002	Short Amplatz Left	
SAL0.75 SH	100cm	HW60AL072P10002	HW70AL072P10002	HW80AL072P10002		
SAL1.0	100cm	HW60AL100P10002	HW70AL100P10002	HW80AL100P10002		
SAL1.0 SH	100cm	HW60AL102P10002	HW70AL102P10002	HW80AL102P10002		
SAL1.5	100cm	HW60AL150P10002	HW70AL150P10002	HW80AL150P10002		
SAL1.5 SH	100cm	HW60AL152P10002	HW70AL152P10002	HW80AL152P10002		
AR1.0	100cm	HW60AR100P10000	HW70AR100P10000	HW80AR100P10000	Amplatz Right	
AR1.0 SH	100cm	HW60AR102P10000	HW70AR102P10000	HW80AR102P10000		
AR2.0	100cm	HW60AR200P10000	HW70AR200P10000	HW80AR200P10000		
AR2.0 SH	100cm	HW60AR202P10000	HW70AR202P10000	HW80AR202P10000		
PB3.0	100cm	HW60PB300P10000	HW70PB300P10000	HW80PB300P10000	Power Backup	
PB3.0 SH	100cm	HW60PB302P10000	HW70PB302P10000	HW80PB302P10000		
PB3.5	100cm	HW60PB350P10000	HW70PB350P10000	HW80PB350P10000		
PB3.5 SH	100cm	HW60PB352P10000	HW70PB352P10000	HW80PB352P10000		
PB4.0	100cm	HW60PB400P10000	HW70PB400P10000	HW80PB400P10000		
PB4.0 SH	100cm	HW60PB402P10000	HW70PB402P10000	HW80PB402P10000		
SPB3.0	100cm	HW60SP300P10000	HW70SP300P10000	HW80SP300P10000	Super Power Backup	
SPB3.0 SH	100cm	HW60SP302P10000	HW70SP302P10000	HW80SP302P10000		
SPB3.5	100cm	HW60SP350P10000	HW70SP350P10000	HW80SP350P10000		
SPB3.5 SH	100cm	HW60SP352P10000	HW70SP352P10000	HW80SP352P10000		
SPB3.75	100cm	HW60SP370P10000	HW70SP370P10000	HW80SP370P10000		
SPB3.75 SH	100cm	HW60SP372P10000	HW70SP372P10000	HW80SP372P10000		
SPB4.0	100cm	HW60SP400P10000	HW70SP400P10000	HW80SP400P10000		
SPB4.0 SH	100cm	HW60SP402P10000	HW70SP402P10000	HW80SP402P10000		
RB0.75	100cm	HW60RB070P10000	HW70RB070P10000	HW80RB070P10000	RCA Backup	
RB0.75 SH	100cm	HW60RB072P10000	HW70RB072P10000	HW80RB072P10000		
RB1.0	100cm	HW60RB100P10000	HW70RB100P10000	HW80RB100P10000		
RB1.0 SH	100cm	HW60RB102P10000	HW70RB102P10000	HW80RB102P10000		
SC3.5	100cm	HW60SC350P10000	HW70SC350P10000	HW80SC350P10000	Special Curve	
SC3.5 SH	100cm	HW60SC352P10000	HW70SC352P10000	HW80SC352P10000		
SC4.0	100cm	HW60SC400P10000	HW70SC400P10000	HW80SC400P10000		
SC4.0 SH	100cm	HW60SC402P10000	HW70SC402P10000	HW80SC402P10000		
HS01	100cm	HW60HS010P10000	HW70HS010P10000	HW80HS010P10000	Hockey Stick	
HS01 SH	100cm	HW60HS012P10000	HW70HS012P10000	HW80HS012P10000		
HS02	100cm	HW60HS020P10000	HW70HS020P10000	HW80HS020P10000		
HS02 SH	100cm	HW60HS022P10000	HW70HS022P10000	HW80HS022P10000		
IM	100cm	HW60IM000P10000	HW70IM000P10000	HW80IM000P10000	Internal Mammary	
IM SH	100cm	HW60IM002P10000	HW70IM002P10000	HW80IM002P10000		
MP01	100cm	HW60MP010P10000	HW70MP010P10000	HW80MP010P10000	Multi Purpose	
MP01 SH	100cm	HW60MP012P10000	HW70MP012P10000	HW80MP012P10000		
MP02	100cm	HW60MP020P10000	HW70MP020P10000	HW80MP020P10000		
MP02 SH	100cm	HW60MP022P10000	HW70MP022P10000	HW80MP022P10000		



ASAHI Hyperion

GUIDING CATHETER

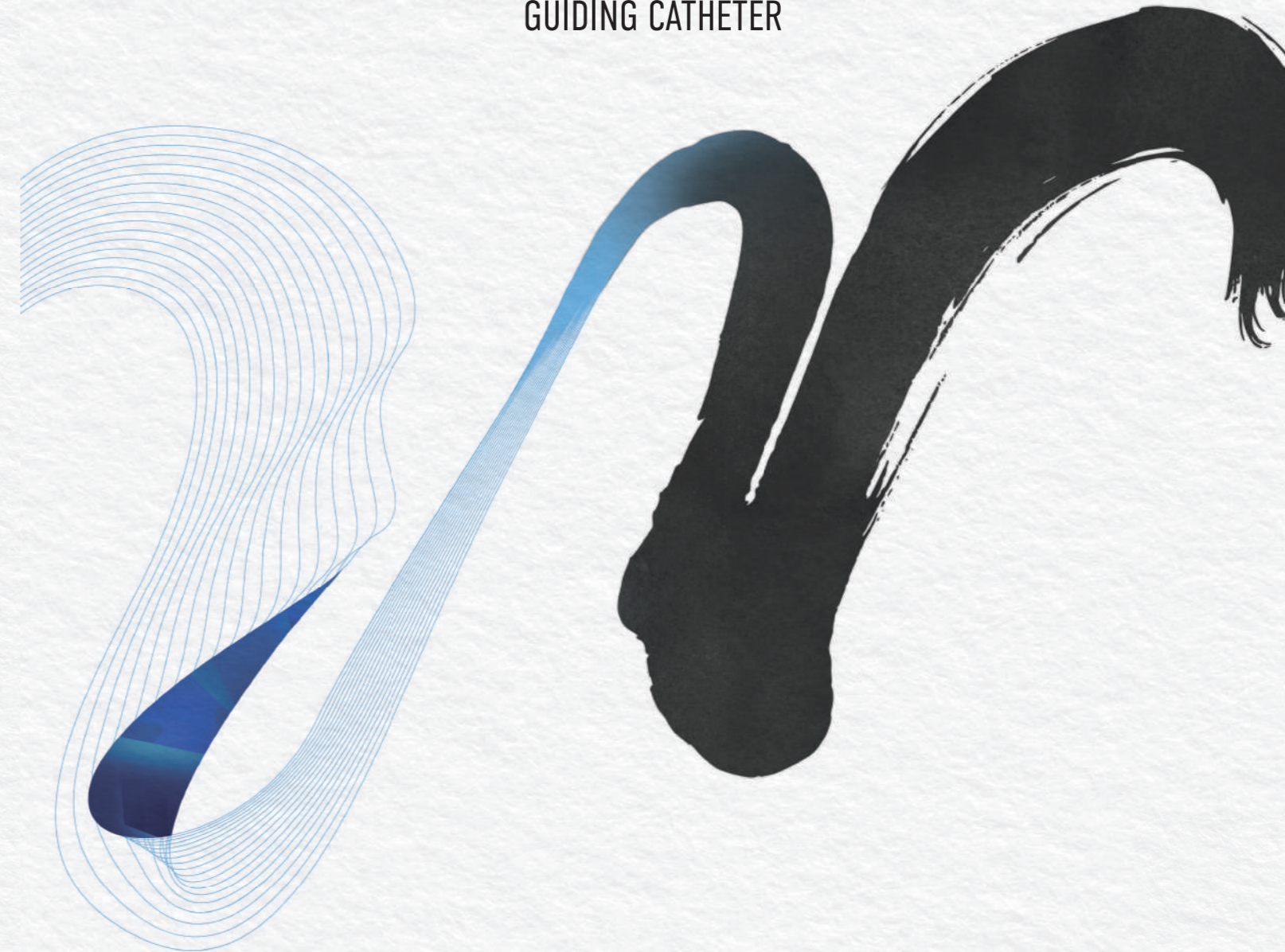

Not all the shapes and sizes are available in all markets. Please check availability with your sales representative.

Some materials/graphics were designed and produced by Vascular Perspective Ltd. www.vascularperspectives.com.

© 2022 ASAHI INTECC CO., LTD. "ASAHI" and "Hyperion" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries. Coated with SLIP-COAT® coating

Manufactured by:

Contact:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
Contact phone number: +81-561-86-9101

www.asahi-intecc.com



ASAHI product information


P22020_B_NOCE_1.0

YOUR DREAMS. WOVEN TOGETHER.

ASAHI Hyperion

Creating a stable environment for your PCI procedure

The universal guide catheter that can navigate difficult PCI conditions has arrived. *The solution lies in here.*

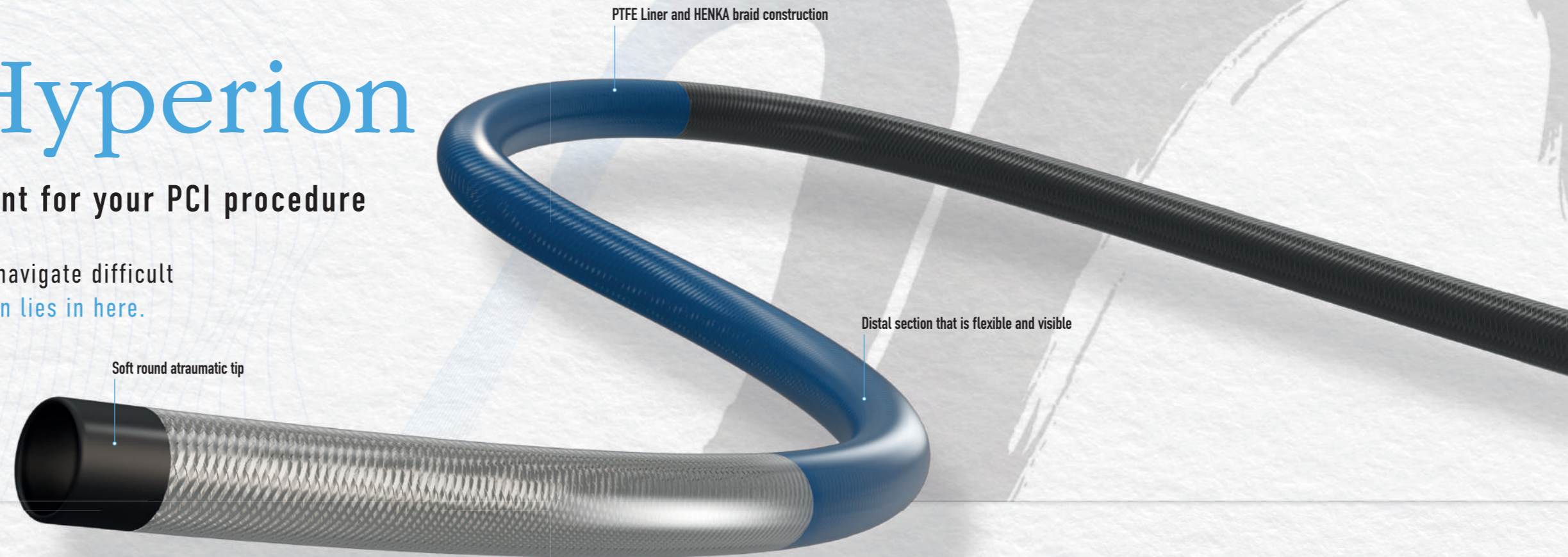


Inner Lumen

The reasonable specifications for your stable PCI procedure.

Enough room for your needs

The large inner diameter diminishes the worry about friction between devices especially with 6F KBT, and provides more effective contrast visualization.

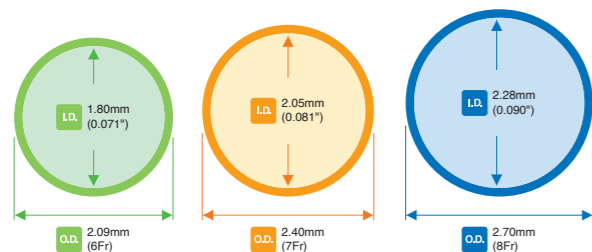

PTFE liner


Smooth inner lumen reduces device insertion resistance.

Flexible Tip

Safe engagement and correct positioning made possible with our flexible-tip technology.


Material

Safe engagement is made possible due to the Urethane tip, the same material used for the tip of ASAHI Corsair, providing both flexibility and high visualization.

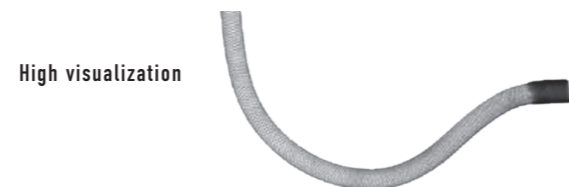

All round processed tip

Processed round tip walls result in atraumatic vessel engagement.

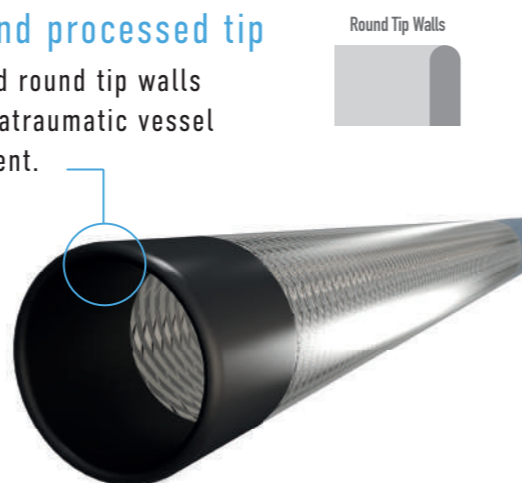

PTFE Liner and HENKA braid construction

Distal section that is flexible and visible

Soft round atraumatic tip

Backup

Long-lasting backup even during challenging procedures.

Hyper Shaft

ASAHI INTECC's unique processing technologies prevent reduction of maneuverability and backup due to heat and blood moisture during the procedure.


Reinforced 2nd curve

A reinforced 2nd curve and excellent shape retention are key factors for high backup.

ASAHI INTECC-original 3rd curve

ASAHI Hyperion's unique 3rd curve at the proximal part of the 2nd curve, improves backup by increasing the contact area between the catheter and the aortic wall.

*For JL, JLST, AL and SAL.

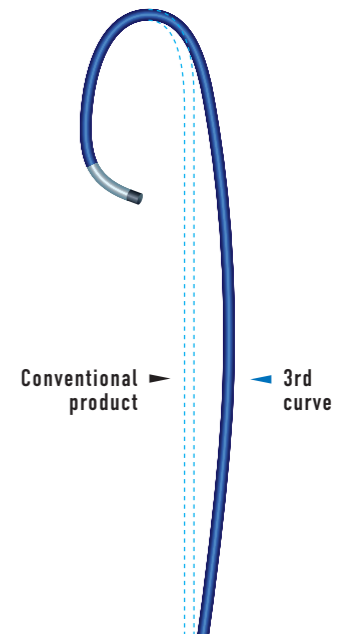


Maneuverability

Provides superior maneuverability in tortuous vessels and prevents the catheter from jumping out of the ostium when engaging.

HENKA-Braid

ASAHI INTECC's unique braiding technology provides well-balanced shaft with reinforced supportive shaft that gradually becomes softer towards the tip.

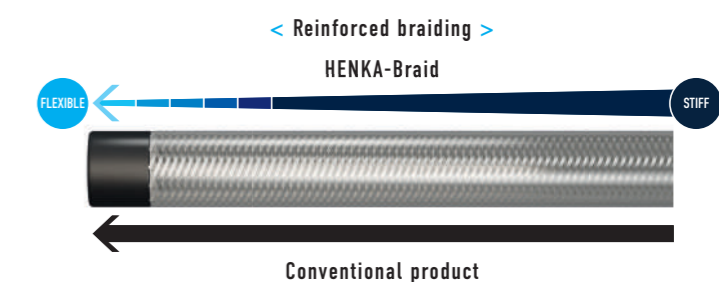

*HENKA is a Japanese word for transformation