

Occlutech Delivery Set

Product specifications

ODS

Perfecting Performance



Occlutech Delivery Set



45° Length 800 mm

Article No.	Size [F]	Length Sheath [mm]	Length Loader [mm]	Dilator Tip [F]	Hemo-Valve [F]	Distal Angle
51DS007	7	800	120	25	70	45°
51DS008	8	800	160	25	70	45°
51DS009	9	800	160	25	70	45°
51DS010	10	800	160	25	70	45°
51DS011	11	800	200	25	70	45°
51DS012	12	800	200	25	70	45°
51DS014	14	800	200	25	70	45°



45° Length 1100 mm

Article No.	Size [F]	Length Sheath [mm]	Length Loader [mm]	Dilator Tip [F]	Hemo-Valve [F]	Distal Angle
51ES006	6	1100	120	25	70	45°
51ES007	7	1100	120	25	70	45°
51ES008	8	1100	160	25	70	45°
51ES009	9	1100	160	25	70	45°
51ES010	10	1100	160	25	70	45°
51ES011	11	1100	200	25	70	45°



180° Length 800 mm

Article No.	Size [F]	Length Sheath [mm]	Length Loader [mm]	Dilator Tip [F]	Hemo-Valve [F]	Distal Angle
51US006	6	800	120	25	70	180°
51US007	7	800	120	25	70	180°
51US008	8	800	160	25	70	180°
51US009	9	800	160	25	70	180°



ODS



Occlutech International AB
Landskronavägen 2, SE-252 32 Helsingborg, Sweden
order@occlutech.com, www.occlutech.com



Perfect Delivery

Features & Benefits

A. Dilator

The dilator (A) is made from smooth conformable low friction radio-opaque material. The dilator is inserted into the sheath (B) to form a supporting inner body. The inner lumen of the dilator has 0.035" guide wire compatibility.

B. Delivery sheath



The sheath (B) consists of a smooth outer layer and a low friction inner layer reinforced with metal braiding. The braid improves the strength of the sheath and provides radiopacity. The sheath has a proximal female connector for connection to the dilator (A) or to the loader (C).

C. Transparent loader

The loader (C) is made of low friction material and is transparent to allow inspection for air bubbles after pulling the occluder into the loader.

D. Y-connector with hemostasis valve

The Y-connector (D) is a commercially available connector with a side-port and hemostasis valve. The valve allows for leak-proof passage of the Occluder delivery wire. The side-port permits flushing before and during implantation.

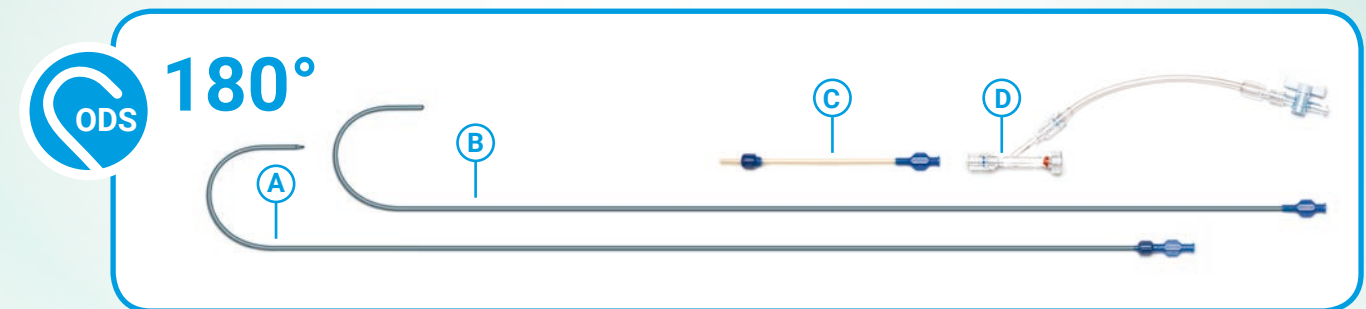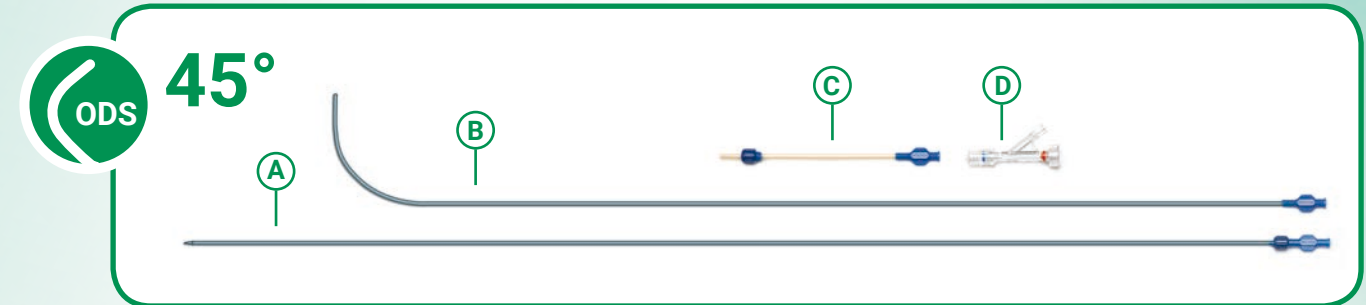

Occlutech Delivery Set

